

OLYMPUS

Your Vision, Our Future

Evolutiontip

CYF-5 FLEXIBLE CYSTOSCOPE



An evolution in smoothness with a refined tip design.



SMOOTH INSERTION WITH EVOLUTIONTIP


Evolutiontip for smooth insertion

The Evolutiontip is a tapered tip that allows for easier insertion. This new design enables smooth entry without sacrificing image quality or treatment capabilities, making the CYF-5 an ideal scope for transurethral procedures.

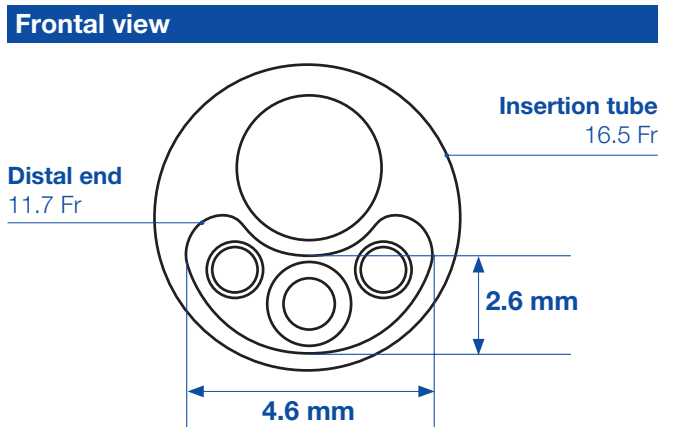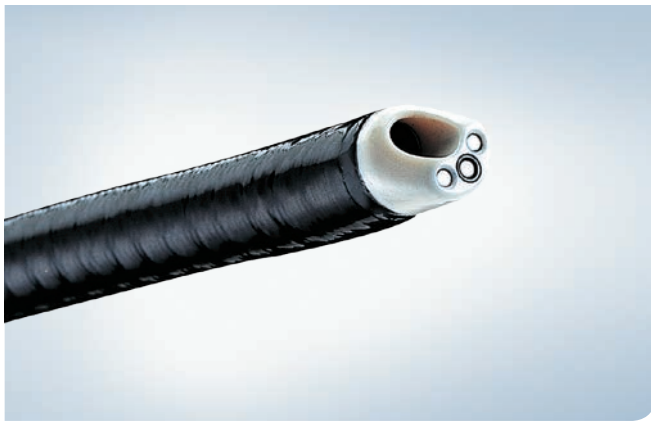

Evolutiontip and matching angulation direction

By matching the angulation direction with the tapered scope end, insertion is facilitated more smoothly.

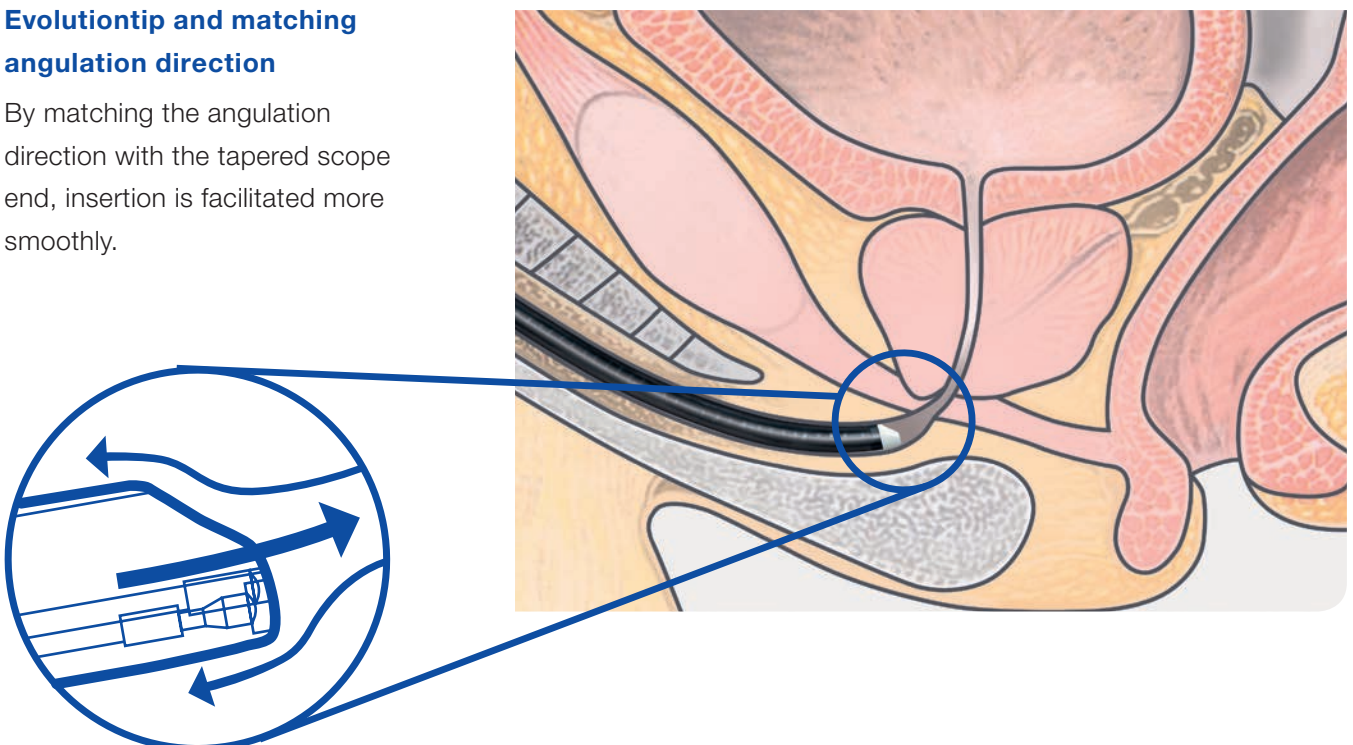

Excellent observation capabilities

- **Improved image quality.**
Newly designed optical system allows for better resolution, especially in the peripheral view.
- **Up 210°/down 120° angulation.** Easier observation and treatment around the ostium and bladder neck.
- **120° field of view.** Both ostia can be observed simultaneously.



Miniature light source

The rechargeable miniature light source with LED (EndoLED WA91502A) transforms the CYF-5 into a portable fibroscope. Perform high-quality examinations anywhere within your medical facility.




Ergonomic design and aspiration function

- **2.4 mm channel.** The large 2.4 mm (7.2 Fr) channel allows effective perfusion with 6 Fr instruments.
- **Aspiration.** One-handed control of the aspiration function makes it easy to keep the field of view clear.
- **Ergonomic design.** Angulation can be locked, allowing single-handed control of the instrument while maintaining the same view.



CYF-5/CYF-5A/CYF-5R

System chart

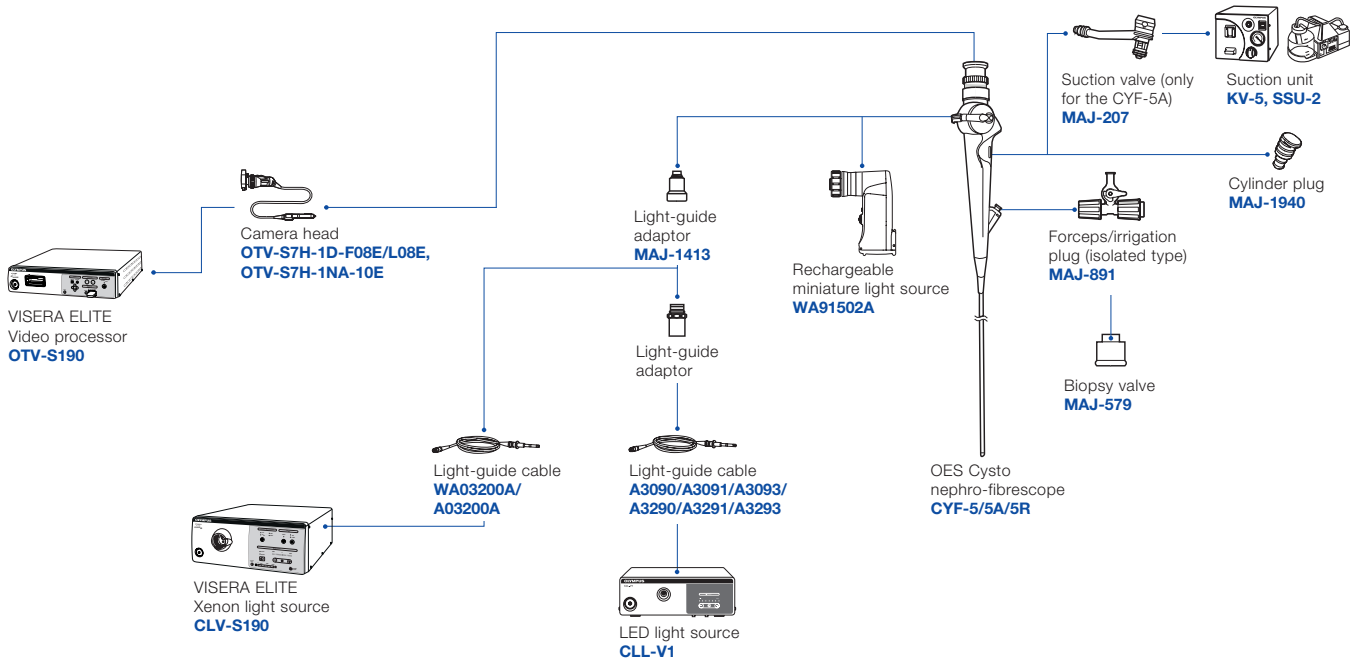

Standard set

WA20810A	CYF-5	Flexible cysto-nephroscope	1	A028958	MAJ-207	Semi-disposable suction valve (CYF-5A only)	10
WA20811A	CYF-5A	Flexible cysto-nephroscope with aspiration	1	026799	BW-15B	Channel brush	1
N4473230	CYF-5R	Flexible cysto-nephroscope with reverse angulation	1	026441	MH-507	Channel-opening cleaning brush	1
N1773930	MAJ-1413	Light-guide adaptor	1	027706	MB-156	ETO cap	1
A02902A	FG-53SX-1	Grasping forceps	1	N1029300	MAJ-1077	Suction cleaning adaptor (CYF-5A only)	1
A02917A	MAJ-891	Forceps/irrigation plug (isolation type)	1				

Specifications

Optical system	Field of view	120°
	Direction of view	Forward-viewing
	Depth of field	3–50 mm
Insertion section	Outer diameter of distal end	11.7 Fr (4.6 mm × 2.6 mm)
	Outer diameter of insertion tube	16.5 Fr (5.5 mm)
Instrument channel	Diameter of inner channel	7.2 Fr (2.4 mm)
	Compatible Olympus EndoTherapy accessories	EndoTherapy accessories for 6 Fr (2 mm) channel
	Minimal visible distance	5 mm
	Bending section	Angulation range
Total length		700 mm

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS

OLYMPUS EUROPA SE & CO. KG

Postbox 10 49 08, 20034 Hamburg, Germany
Wendenstrasse 14–18, 20097 Hamburg, Germany
Phone: +49 (0)40 237 730, Fax: +49 (0)40 230 761
www.olympus-europa.com